

SEQUENTIAL CIRCULATOR 2004

HCPCS CODE E0651

A Gradient, Sequential, Pneumatic Compression Device for...

LYMPHEDEMA VENOUS INSUFFICIENCY WOUND CARE

INTENDED USE


The **Model SC-2004 Sequential Circulator** is a gradient, sequential, pneumatic compression device, intended for the primary or adjunctive treatment of primary or secondary Lymphedema. The device is additionally intended for alternate treatment of chronic venous stasis ulcers and associated venous insufficiency, as well as general treatment of swelling of the extremities. The device is intended for home or hospital use.

*Latest Clinical Studies using
BIO COMPRESSION SYSTEMS'
SEQUENTIAL CIRCULATORS
are available upon request*



Bio Compression Systems Sequential Circulator Model 2004...

- *Reduces Pain*
- *Mimics the Lymphatic System*
- *Promotes Lymphatic Flow by Moving Fluid in the Proper Physiological Direction (distal to proximal)*
- *Stimulates Circulation to Promote Wound Healing*
- *Delivers Continuous Flow*
- *Bilateral*
- *Provides Effective Prophylaxis for Venous Thrombosis*
- *Lightweight • Portable • Easy to Use • Cost-Effective*

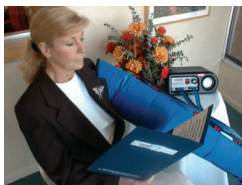

Arm and Hand Therapy

Bio Compression Systems' Biocomfort Garments ...

Made of high quality 200 Denier Nylon Oxford coated with 3 mils of Polyurethane

Available in a variety of sizes and styles including

OVER THE SHOULDER • ADJUSTABLE • CUSTOM • AMPUTEE



Foot and Leg Therapy

BIOCOMFORT GARMENT SIZE CHART ON BACK

CONTRAINDICATIONS

Pneumatic compression is contraindicated for the patients with the following: 1) Congestive heart failure; 2) Deep Vein Thrombosis; 3) Inflammatory phlebitis or episodes of pulmonary embolism; 4) Infections in the limb, including cellulitis, without appropriate antibiotic coverage; 5) Presence of lymphangiosarcoma
WARNING: Federal law restricts this device to sale by or on the order of a physician.



#1 in Compression Therapy

120 West Commercial Avenue, Moonachie, NJ 07074

Phone: 800-888-0908 Fax: 201-939-4503

E-mail: biosystems@biocompression.com

Website: www.biocompression.com

SPECIFICATIONS

MODEL 2004

HCPSC CODE E0651

Segmental, pneumatic compression providing pre-determined gradient pressure.

Dimensions H x W x D: 5.5 x 12 x 8
 Weight: 6 lbs.
 Pressure Range: 25-100mm Hg
 Electrical Rating: 120 VAC, 60 HZ, .5 A
 Cycle Time: 18 secs per chamber
 Inflation : 72 secs.
 Deflation: 18 secs.
 Warranty: 3 years
 Bi-lateral operation

BIOCOMFORT GARMENTS SIZE CHART *Now in* INCHES/CENTIMETERS

<u>Upper Extremity</u>	<i>1</i>	<i>2</i>	<i>3</i>	<i>4</i>	<i>5</i>	<i>6</i>	<i>6</i> Axilla to Fingertips Length
	<u>Axilla</u>	<u>Bicep</u>	<u>Elbow</u>	<u>Forearm</u>	<u>Wrist</u>		
<u><i>Small Arm</i></u>							
GS-3035-S	27"/69.6	26"/66.0	24"/61.0	23"/58.4	19.5"/49.5		26.5"/67.3
<u><i>Medium Arm</i></u>							
GS-3035-M	27"/69.6	26"/66.0	24"/61.0	23"/58.4	19.5"/49.5		31.0"/78.7
<u><i>Large Arm</i></u>							
GS-3035-L	27"/69.6	26"/66.0	24"/61.0	23"/58.4	19.5"/49.5		35.5"/91.2
<u><i>Arm w/ Shldr</i></u>							
GS-3035 SH	27"/69.6	26"/66.0	24"/61.0	23"/58.4	19.5"/49.5		31+8"/78.7+20.3
<u>Lower Extremity</u>	<i>1</i>	<i>2</i>	<i>3</i>	<i>4</i>	<i>5</i>	<i>6</i>	<i>7</i>
	<u>Thigh below Crotch</u>	<u>Mid Thigh</u>	<u>Mid Knee</u>	<u>Mid Calf</u>	<u>Ankle</u>	<u>Arch (Instep)</u>	<u>Crotch to Heel</u>
<u><i>Half Leg</i></u>	—	—	29"/73.7	25.5"/64.8	24.5"/62.2	21.5"/54.6	—
GS-3045-H							
<u><i>Small Leg</i></u>							
GS-3045-S	32"/81.3	31"/78.7	25"/63.5	24.5"/62.2	24.5"/62.2	21.5"/54.6	27"/69.6
<u><i>Medium Leg</i></u>							
GS-3045-M	32"/81.3	31"/78.7	27"/68.6	25"/63.5	24.5"/62.2	21.5"/54.6	32"/81.3
<u><i>Large Leg</i></u>							
GS-3045-L	32"/81.3	31"/78.7	29"/73.7	25.5/ 64.8	24.5"/62.2	21.5"/ 54.6	36"/91.4
<u><i>Half Leg Adjustable (Expanded Position)</i></u>	—	—	37.5"/95.3	32"/81.3	24.5/62.2	21.5"/54.6	—
GA-3045-H							
<u><i>Small Leg Adjustable (Expanded Position)</i></u>	45.5"/115.6	42"/106.7	33.5"/85.1	31"/78.8	24.5"/62.2	21.5"/54.6	27"/69.6
GA-3045-S							
<u><i>Medium Leg Adjustable (Expanded Position)</i></u>	45.5"/115.6	42"/106.7	35.5"/90.2	31.5"/80	24.5"/62.2	21.5"/54.6	32"/81.3
GA-3045-M							
<u><i>Large Leg Adjustable (Expanded Position)</i></u>	45.5"/115.6	42"/106.7	37.5"/95.3	32"/81.3	24.5"/62.2	21.5"/54.6	36"/91.4
GA-3045-L							